

PRODUCT DATA SHEET

Published: 2026-07-30

LAUDA Variocool VC 1200

Process thermostat 230 V; 50 Hz

Part Number: L000923

Features

- Process thermostat suitable for use with non-flammable heat transfer liquids
- Coloured TFT display for simultaneous indication of actual & set values and graphic illustration of the temperature profile
- Clear text menu navigation, six selectable languages DE, EN, FR, ES, IT, RU
- Easy input via cursor and soft keys
- Fully electronic continuous controller with PID action
- Electronic level indication and low level alarm
- Powerful pressure pump
- USB interface as standard
- Remote fault indication through floating contact
- Upgradeable with an interface module (analogue module, contact module, RS 232/485 module, Profibus, Ethernet-USB module)
- Integrated programmer with max. 150 segments, splitable in 5 programmes
- Adjustable bypass for pressure limiting
- Filler opening on top, drain tap on the backside
- SmartCool system for energy-saving digital cooling management including compressor on-off control
- Utilises traditional refrigerants (HFCs) in accordance with European legislation to control F-gases (EU) 573/2024
- Condenser cooling Air
- Operates with non flammable liquids (water, water/glycol)



Working temperature min.
-20 °C



Working temperature max.
80 °C

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser

PRODUCT DATA SHEET

Published: 2026-07-30


LAUDA Variocool VC 1200

Process thermostat 230 V; 50 Hz

Part Number: L000923

Technical Features (according to DIN 12876)

Working temperature range	-20 ... 80 °C
Ambient temperature range	5 ... 40 °C
Temperature stability	0.05 ± K
Heater power max.	2.25 kW
Power consumption max.	3.3 kW
Current max.	14 A
Pump Pressure max.	3,2 bar
Pump flow rate max. (pressure)	37 L/min
In / Outlet connection thread (outside)	G 3/4"
Pressure adjustment	bypass
Filling volume min.	8 L
Filling volume max.	15 L
Overall dimensions (WxDxH)	450 x 550 x 790 mm
Weight	51 kg
Refrigerant stage 1	R-449A (GWP 1397); 0.500 kg; 0.7 t CO2-eq
Power supply	230 V; 50 Hz
Power plug	Power cord with angled plug (CEE7/7)

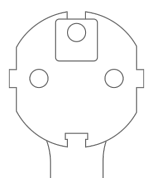

Power cord with angled plug (CEE7/7)

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser

PRODUCT DATA SHEET

Published: 2026-07-30

LAUDA Variocool VC 1200

Process thermostat 230 V; 50 Hz

Part Number: L000923

Temperature	Heat transfer liquid	Cooling Capacity 50 Hz
20 °C	Ethanol	1.12 kW
10 °C	Ethanol	0.92 kW
0 °C	Ethanol	0.62 kW
-10 °C	Ethanol	0.32 kW
-20 °C	Ethanol	0.06 kW

Standard accessories

- 2 nipples 3/4" with screw cap G3/4 for pump connectors

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser

PRODUCT DATA SHEET

Published: 2026-07-30

LAUDA Variocool VC 1200

Process thermostat 230 V; 50 Hz

Part Number: L000923

PUMP CHARACTERISTIC Heat transfer liquid: Water

Pressure bar



- 1 0.9 bar, 28 L/min
- 2 3.2 bar, 37 L/min
- 3 4.8 bar, 37 L/min
- 4 5.0 bar, 60 L/min

Flow rate L/min

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser